

(19)



Europäisches Patentamt

European Patent Office

Office européen des brevets



(11)

EP 1 050 587 A2

(12)

EUROPEAN PATENT APPLICATION

(43) Date of publication:
08.11.2000 Bulletin 2000/45

(51) Int. Cl.⁷: **C12Q 1/68**

(21) Application number: **00109078.6**

(22) Date of filing: **02.05.2000**

(84) Designated Contracting States:
**AT BE CH CY DE DK ES FI FR GB GR IE IT LI LU
MC NL PT SE**
Designated Extension States:
AL LT LV MK RO SI

(30) Priority: **03.05.1999 US 304452**

(71) Applicant: **QIAGEN GmbH
40724 Hilden (DE)**

(72) Inventors:
• **Missel, Andreas
40215 Düsseldorf (DE)**

• **Loeffert, Dirk
40589 Düsseldorf (DE)**
• **Kang, Jie
40822 Mettmann (DE)**
• **Korfhage, Christian
40674 Langenfeld (DE)**

(74) Representative:
**Mahlbacher, Volker et al
QIAGEN GmbH,
Intellectual Property & Licensing,
Max-Volmer-Strasse 4
40724 Hilden (DE)**

(54) Generation and amplification of nucleic acids from ribonucleic acids

(57) Novel compositions and methods useful for the generation of nucleic acids from a ribonucleic acid template and further nucleic acid replication are disclosed. It is shown that the generation and amplification of nucleic acids by methods that utilize two or more different polymerases, such as reverse transcriptase-polymerase chain reaction (RT-PCR), are dramatically more sensitive and efficient in the presence of a homopolymeric nucleic acid. Homopolymeric nucleic acids have been found to reduce or negate the inhibitory effect reverse transcriptases have on DNA polymerase activity. It is demonstrated that this inhibition-relieving effect of homopolymeric nucleic acids is general in nature; independent of the chemical species of homopolymer used, or the chemical composition of the polymerization reaction mixture.

EP 1 050 587 A2

Description

Field of the Invention

5 [0001] The present invention is in the fields of molecular biology. The present invention is directed to novel compositions and methods useful for the generation of nucleic acids from a ribonucleic acid template and further nucleic acid replication. Specifically, the invention is directed to the generation and amplification of nucleic acids by reverse transcriptase-polymerase chain reaction (RT-PCR).

10 Background of the Invention

[0002] The detection, analysis, transcription, and amplification of nucleic acids are the most important procedures in modern molecular biology. The application of such procedures for RNA analysis is especially important in the investigation of gene expression, diagnosis of infectious agents or genetic diseases, the generation of cDNA, and analysis of retroviruses, to name but a few applications. The reverse transcription of RNA, followed by polymerase chain reaction amplification, commonly referred to as RT-PCR or RNA-PCR, has become widely used for the detection and quantification of RNA.

[0003] The RT-PCR procedure involves two separate molecular syntheses: (i) the synthesis of cDNA from an RNA template; and (ii) the replication of the newly synthesized cDNA through PCR amplification. RT-PCR may be performed under three general protocols: (1) uncoupled RT-PCR, also referred to as two-step RT-PCR; (2) single enzyme coupled RT-PCR (coupled RT-PCR is also referred to as one-step RT-PCR or continuous RT-PCR), in which a single polymerase is used for both the cDNA generation from RNA as well as subsequent DNA amplification; and (3) two (or more) enzyme coupled RT-PCR, in which at least two separate polymerases are used for initial cDNA synthesis and subsequent replication.

25 [0004] In uncoupled RT-PCR, reverse transcription is performed as an independent step using buffer and reaction conditions optimal for reverse transcriptase activity. Following cDNA synthesis, an aliquot of the RT reaction product is used as template for PCR amplification with a thermostable DNA Polymerase, such as *Taq* DNA Polymerase, under conditions optimal for PCR amplification.

[0005] In coupled RT-PCR, reverse transcription and PCR amplification are combined into a single reaction mixture. Single enzyme RT-PCR utilizes the reverse transcriptase activity of some DNA polymerases, such as *Taq* DNA Polymerase and *Tth* DNA polymerase, whereas two-enzyme RT-PCR typically uses a retroviral or bacterial reverse transcriptase (e. g. AMV-RT, MMLV-RT, HIV-RT, EIAV-RT, RAV2-RT, *Carboxydotherrnus hydrogeniformans* DNA Polymerase or a mutant, variant or derivative thereof), and a thermostable DNA polymerase (e. g. *Taq*, *Tbr*, *Tth*, *Tih*, *Tfi*, *Tfl*, *Pfu*, *Pwo*, *Kod*, *VENT*, *DEEPVENT*, *Tma*, *Tne*, *Bst*, *Pho*, *Sac*, *Sso*, *ES4* and others or a mutant, variant or derivative thereof).

[0006] Coupled RT-PCR provides numerous advantages over uncoupled RT-PCR. Coupled RT-PCR requires less handling of the reaction mixture reagents and nucleic acid products than uncoupled RT-PCR (e.g., opening of the reaction tube for component or enzyme addition in between the two reaction steps), and is therefore less labor intensive, reducing the required number of person hours. Coupled RT-PCR also requires less sample, and reduces the risk of contamination (Sellner and Turbett, 1998).

[0007] Single enzyme coupled RT-PCR, is the simplest RT-PCR procedure to date. This system is expensive to perform, however, due to the amount of DNA polymerase required. In addition, the single enzyme coupled RT-PCR method has been found to be less sensitive than uncoupled RT-PCR (Cusi et al., 1994), and limited to polymerizing nucleic acids of less than one kilobase pair (>1 kpb) in length. Two enzyme RT-PCR systems show increased sensitivity over the single enzyme system generally, even when coupled in a single reaction mixture. This effect has been attributed to the higher efficiency of reverse transcriptase in comparison to the reverse transcriptase activity of DNA polymerases (Sellner and Turbett, 1998).

[0008] Although the two enzyme coupled RT-PCR system is more sensitive than the uncoupled protocol, reverse transcriptase has been found to interfere directly with DNA polymerase during the replication of the cDNA, thus reducing the sensitivity and efficiency of this technique (Sellner et al., 1992; Aatsinki et al., 1994; Mallet et al., 1995). A variety of solutions to overcome the inhibitory activity of reverse transcriptase on DNA polymerase have been tried, including: increasing the amount of template RNA, increasing the ratio of DNA polymerase to reverse transcriptase, adding modifier reagents that may reduce the inhibitory effect of reverse transcriptase on DNA polymerase (e.g., non-homologous tRNA, T4 gene 32 protein, sulfur or acetate-containing molecules), and heat-inactivation of the reverse transcriptase before the addition of DNA polymerase.

55 [0009] All of these modified RT-PCR methods have significant drawbacks, however. Increasing the amount of template RNA is not possible in cases where only limited amounts of sample are available. Individual optimization of the ratio of reverse transcriptase to DNA polymerase is not practicable for ready-to-use reagent kits for one-step RT-PCR.

The net effect of currently proposed modifier reagents to relieve reverse transcriptase inhibition of DNA polymerization is controversial and in dispute: positive effects due to these reagents are highly dependent on RNA template amounts, RNA composition, or may require specific reverse transcriptase-DNA polymerase combinations (see, for example, Chandler et al., 1998). Finally, heat inactivation of the reverse transcriptase before the addition of the DNA polymerase negates the advantages of the coupled RT-PCR and carries with it all the disadvantages of uncoupled RT-PCR systems discussed earlier.

[0010] Because of the importance of RT-PCR applications, a coupled RT-PCR system, in the form of a generalized ready-to-use composition, which exhibits high sensitivity, requires a small amount of initial sample, reduces the amount of practitioner manipulation, minimizes the risks of contamination, minimizes the expense of reagents, is not restricted to the use of specific reaction buffers, and maximizes the amount of nucleic acid end product is needed in the art.

Summary of the Invention

[0011] The present invention is directed to novel compositions and methods useful for the generation of nucleic acids from a ribonucleic acid template and further nucleic acid replication. Specifically, the invention is directed to the generation and amplification of nucleic acids by reverse transcriptase-polymerase chain reaction (RT-PCR), utilizing two or more different polymerases in the presence of a homopolymeric nucleic acid. The presence of the homopolymeric nucleic acid serves to negate inhibition of the DNA polymerase, which occurs when a composition comprising two or more polymerases, at least one of which having reverse transcriptase activity, is used. The compositions and methods of present invention are used to generate and replicate DNA molecules complementary to a portion of an RNA template.

[0012] The present invention is directed generally to all reaction mixtures that can be used in the generation and replication of a nucleic acid from an RNA template. One embodiment of the present invention is directed to a composition comprising a homopolymeric nucleic acid and one or more of the following reaction components: a reverse transcriptase, a DNA polymerase, one or more oligonucleotide primers, any one or more nucleotide base, an appropriate buffering agent, a salt solution, or other additives useful in RT-PCR.

[0013] The present invention offers several advantages compared to known methods for generating cDNA from an RNA target, including, but not limited to:

- permitting coupled RT-PCR, involving one reaction solution and reduced handling of reagents and products;
- permitting use of a small initial sample of RNA template;
- permitting use of a wide range of different reverse transcriptases;
- permitting use of a wide range of different DNA polymerases;
- permitting use of a wide range of different reaction mixture buffer and salt solutions, including buffers containing additional specialized reaction additives (e.g., sulfate- and acetate-containing compounds);
- reducing adverse effects on the specificity and product yield (observed with tRNA additives);
- operating with reagents that are commercially available, easily synthesized, and inexpensive.

The invention thus facilitates the rapid and efficient generation of nucleic acid molecules from a sample containing ribonucleic acids (RNA) as well as the detection and quantitation of RNA molecules. The invention also is useful in the rapid production and amplification of cDNAs that may be used for a variety of industrial, medical and forensic purposes.

[0014] The homopolymeric nucleic acid can be a ribonucleic acid homopolymer such as poly(A), poly(C), poly(G) or poly(U), or a deoxyribonucleic acid homopolymer such as poly(dA), poly(dC), poly(dG), poly(dT) or poly(dI). Preferred homopolymers are poly(A), poly(C), poly(G), poly(U), poly(dA), poly(dC), poly(dG), and poly(dI), or a mixture of two or more of these homopolymers. The homopolymer used for RT-PCR is an oligonucleotide of two or more bases in length, preferably between 100 to 5000 bases in length, and more preferably between 200 to 2000 bases in length. The amount of nucleic acid homopolymer used for RT-PCR may be greater than 10 ng, preferably between 10 ng to 2000 ng, more preferably between 100 ng to 1000 ng, and most preferably between 200 ng to 600 ng. Preferably, the final concentration of nucleic acid homopolymer in the RT-PCR is between 200 ng/ml to 40000 ng/ml, more preferably between 2000 ng/ml to 20000 ng/ml, and most preferably between 4000 ng/ml to 12000 ng/ml. The present invention includes compositions containing nucleic acid homopolymer in higher concentrations as the preferred final concentrations given above, but which are intended to achieve the preferred final concentrations by dilution, e.g. by combination with further components of the final reaction mixture.

[0015] The reverse transcriptase may be any polymerase that exhibits reverse transcriptase activity. Preferably, the reverse transcriptase is AMV-RT, MMLV-RT, HIV-RT, EIAV-RT, RAV2-RT, *C. hydrogeniformans* DNA Polymerase, SuperScript I, SuperScript II, or mutants, variants and derivatives thereof having reverse transcriptase activity.

[0016] As used herein, mutants, variants and derivatives refer to all permutations of a chemical species, which may exist or be produced, that still retain the definitive chemical activity of that chemical species: Examples include, but are

not limited to compounds that may be detectably labeled, or otherwise modified, thus altering the compound's chemical or physical characteristics.

[0017] The DNA polymerase may be any polymerase capable of replicating a DNA molecule. Preferably, the DNA polymerase is a thermostable polymerase useful in PCR. More preferably, the DNA polymerase is *Taq*, *Tbr*, *Tth*, *Tfi*, *Tfi*, *Pfu*, *Pwo*, *Kod*, *VENT*, *DEEPVENT*, *Tma*, *Tne*, *Bst*, *Pho*, *Sac*, *Sso*, *Poc*, *Pab*, *ES4* or mutants, variants and derivatives thereof having DNA polymerase activity.

[0018] Oligonucleotide primers may be any oligonucleotide of two or more nucleotides in length. Primers may be random primers, homopolymers, or primers specific to a target RNA template (e.g., a sequence specific primer).

[0019] Additional compositional embodiments comprise a homopolymeric nucleic acid and other reaction mixture components such as one or more nucleotides or derivatives thereof. Preferably the nucleotide is a deoxynucleotide triphosphate: dNTP (e.g., dATP, dCTP, dGTP, dTTP, dITP, dUTP, α -thio-dNTPs, biotin-dUTP, fluorescein-dUTP, digoxigenin-dUTP).

[0020] Buffering agents, salt solutions and other additives of the present invention comprise those solutions useful in RT-PCR. Preferred buffering agents include, e.g., TRIS, TRICINE, BIS-TRICINE, HEPES, MOPS, TES, TAPS, PIPES, CAPS. Preferred salt solutions include, e.g., potassium chloride, potassium acetate, potassium sulfate, ammonium sulfate, ammonium chloride, ammonium acetate, magnesium chloride, magnesium acetate, magnesium sulfate, manganese chloride, manganese acetate, manganese sulfate, sodium chloride, sodium acetate, lithium chloride, and lithium acetate. Preferred additives include, e.g., DMSO, glycerol, formamide, betain, tetramethylammonium chloride, PEG, Tween 20, NP 40, ectoine, polyols, *E. coli* SSB protein, Phage T4 gene 32 protein, BSA.

[0021] Further embodiments of this invention relate to methods for generating nucleic acids from a ribonucleic acid template and further nucleic acid replication. The method comprises: (a) adding an RNA template to a reaction mixture comprising a reverse transcriptase, a nucleic acid polymerase, or derivatives thereof, and a homopolymeric nucleic acid, and (b) incubating the reaction mixture under conditions sufficient to allow polymerization of a nucleic acid molecule complementary to at least a portion of the RNA template. In a preferred embodiment the method includes replication of the DNA molecule complementary to at least a portion of the RNA template. More preferably the method of DNA replication is polymerase chain reaction (PCR). Most preferably the method comprises coupled reverse transcriptase-polymerase chain reaction (RT-PCR).

[0022] Preferably, the nucleic acid homopolymer in such method is poly(A), poly(C), poly(G), poly(U), poly(dA), poly(dC), poly(dG), or poly(dI), or a mixture of two or more of these homopolymers. The homopolymer is preferably an oligonucleotide of two or more bases in length, preferably between 100 to 5000 bases in length, and more preferably between 200 to 2000 bases in length.

[0023] Preferably, the final concentration of nucleic acid homopolymer in the method of the invention is between 200 ng/ml to 40000 ng/ml, more preferably between 2000 ng/ml to 20000 ng/ml, and most preferably between 4000 ng/ml to 12000 ng/ml.

[0024] In more specific embodiments, RT-PCR occurs between about 23°C to about 100°C. Preferably reverse transcription occurs between about 35°C to about 75°C, followed by PCR occurring between about 60°C to about 95°C. Under the most preferred conditions reverse transcription occurs between about 37°C to about 60°C, denaturation occurs at about 94°C, annealization occurs at about 60°C, and polymerization occurs at about 72°C.

Brief Description of the Drawings

[0025]

Fig. 1 is a photograph of an ethidium bromide (EtdBr)-stained 1% agarose gel, which depicts the reduction of reverse transcriptase inhibition of PCR by homopolymeric poly(A). Lane 1 is a nucleic acid ladder, used as a gel reference marker. Lanes 2 and 3 are controls, which lack homopolymeric poly(A). Lanes 4 and 5 show the nucleic acid products of an RT-PCR reaction in the presence of 100 ng of poly(A). Lanes 6 and 7 show the nucleic acid products of an RT-PCR reaction in the presence of 200 ng of poly(A). Lanes 8 and 9 show the nucleic acid products of an RT-PCR reaction in the presence of 400 ng of poly(A).

Fig. 2 is a photograph of an ethidium bromide (EtdBr)-stained 1% agarose gel, which depicts the reduction of reverse transcriptase inhibition of PCR by; **2a** homopolymeric poly(U), and **2b** homopolymeric poly(C). Lane 1 is a nucleic acid ladder, used as a gel reference marker. Lanes 2 and 3 are controls, which lack a homopolymeric nucleic acid. Lanes 4 and 5 show the nucleic acid products of an RT-PCR reaction in the presence of 100 ng of (**2a**) poly(U), or (**2b**) poly (C). Lanes 6 and 7 show the nucleic acid products of an RT-PCR reaction in the presence of 200 ng of (**2a**) poly(U), or (**2b**) poly (C). Lanes 8 and 9 show the nucleic acid products of an RT-PCR reaction in the presence of 400 ng of (**2a**) poly(U), or (**2b**) poly (C).

Fig. 3 is a photograph of an ethidium bromide (EtdBr)-stained 1% agarose gel, which depicts the reduction of reverse transcriptase inhibition of PCR by homopolymeric poly(dA). Lane 1 is a nucleic acid ladder, used as a gel

reference marker. Lanes 2 through 4 are controls, which lack homopolymeric poly(dA). Lanes 5 through 7 show the nucleic acid products of an RT-PCR reaction in the presence of 200 ng of poly(dA). Lanes 8 through 10 show the nucleic acid products of an RT-PCR reaction in the presence of 400 ng of poly(A).

Fig. 4 is a photograph of an ethidium bromide (EtdBr)-stained 1% agarose gel, which depicts the reduction effects of homopolymeric poly(A) on PCR inhibition involving different reverse transcriptases. Lane 1 is a nucleic acid ladder, used as a gel reference marker. Lanes 2 through 4, 9 through 11, and 16 through 18 are controls, which lack homopolymeric poly(A). Lanes 5 through 7 show the nucleic acid products of an RT-PCR reaction using AMV-RT. Lanes 12 through 14 show the nucleic acid products of an RT-PCR reaction using MMLV-RT. Lanes 19 through 21 show the nucleic acid products of an RT-PCR reaction using MMLV-RT RNase H-. (Lanes 8 and 15 are unused as spacers.)

Fig. 5 is a photograph of an ethidium bromide (EtdBr)-stained 1% agarose gel, which depicts the reduction effects of homopolymeric poly(A) on different DNA polymerases and under different buffer and salt conditions. **Fig. 5a** depicts the RT-PCR nucleic acid products using AmpliTaq DNA polymerase in the absence (lanes 2 through 4) of homopolymeric poly(A), and in the presence (lanes 5 through 7) of homopolymeric poly(A). Lane 1 is a nucleic acid ladder, used as a gel reference marker. **Fig. 5b** depicts the RT-PCR nucleic acid products using PFU DNA polymerase in the absence (lanes 2 through 4) of homopolymeric poly(A), and in the presence (lanes 5 through 7) of homopolymeric poly(A). Lanes 1 and 8 are nucleic acid ladders, used as a gel reference marker.

Fig. 6 is a photograph of an ethidium bromide (EtdBr)-stained 1% agarose gel, which depicts the reduction of reverse transcriptase inhibition of PCR by a homopolymeric nucleic acid compared to that of tRNA and rRNA. Lane 1 is a nucleic acid ladder, used as a gel reference marker. Lanes 3 through 5 are controls, which lack any nucleic acid additive to the standard RT-PCR reaction mixture. Lanes 7 through 9 show the nucleic acid products of an RT-PCR reaction in the presence of homopolymeric poly(A). Lanes 11 through 13 show the nucleic acid products of an RT-PCR reaction in the presence of tRNA. Lanes 15 through 17 show the nucleic acid products of an RT-PCR reaction in the presence of rRNA. (Lanes 2, 6, 10, and 14 are unused as spacers.)

Detailed Description of the Invention

[0026] The present invention is directed to compositions and methods for use in the generation of nucleic acids from an RNA template, and particularly for reverse transcriptase-polymerase chain reaction (RT-PCR) production and analysis of nucleic acids. The invention provides compositions that include a homopolymeric nucleic acid in combination with one or more components useful in the generation of nucleic acids from an RNA template. Such components may include a reverse transcriptase, a DNA polymerase, one or more oligonucleotide primers, any one or more nucleotide bases, an appropriate buffering agent, a salt, or other additives useful in RT or RT-PCR. Methods of the invention utilize homopolymeric nucleic acid in combination with other reagents to generate nucleic acid from an RNA template, and preferably also to replicate the newly synthesized nucleic acid. The compositions and methods of the present invention useful to generate, replicate, analyze, quantitate and otherwise manipulate nucleic acid molecules are most useful in coupled or uncoupled RT-PCR procedures.

[0027] RT-PCR is one molecular manipulation used to generate and replicate a nucleic acid derived from an RNA template. RT-PCR is described herein as an exemplary protocol capable of utilizing the compositions and methods of the present invention without limitation. It is understood by one of ordinary skill in the art that the present invention has utility in other processes, which involve a combination of reverse transcriptase and DNA polymerase activity. RT-PCR involves two separate molecular syntheses: (i) the synthesis of cDNA from an RNA template; and (ii) the replication of the newly synthesized cDNA through PCR amplification.

Reverse Transcriptase-Polymerase Chain Reaction (RT-PCR)

[0028] In RT-PCR, the reaction mixture is first incubated (in an appropriate buffering agent) at a temperature sufficient to allow synthesis of a DNA molecule complementary to at least a portion of an RNA template. Components of a reverse transcription reaction mixture typically include an RNA template, from which the complementary DNA (cDNA) is transcribed; a nucleic acid polymerase that exhibits reverse transcriptase activity; and the appropriate nucleotide building blocks needed for nucleic acid synthesis. For the purposes of this invention, cDNA is defined as any DNA molecule whose nucleic acid sequence is complementary to an RNA molecule. An RNA template is defined as any RNA molecule used to provide a nucleic acid sequence from which a cDNA molecule may be synthesized. The synthesis of cDNA from an RNA template is typically accomplished by utilizing a nucleic acid polymerase that exhibits reverse transcriptase activity. For the purposes of this invention, reverse transcriptase activity refers to the ability of an enzyme to polymerize a cDNA molecule from an RNA template, and reverse transcriptase broadly refers to any enzyme possessing reverse transcriptase activity. Reverse transcription typically occurs in a temperature range from about 20°C to about 75°C, preferably from about 35°C to about 70°C.

[0029] After reverse transcription of an RNA template to produce a cDNA molecule, the cDNA is incubated (in an appropriate buffering agent) under conditions sufficient for replication of the cDNA molecule. The reaction mixture may be the same as that of the previous reverse transcription reaction mixture, as employed in coupled (also called continuous, or one-step) RT-PCR, or the reaction mixture may comprise an aliquot of the previous reverse transcription reaction mixture and may be further modified for nucleic acid replication, as in uncoupled (or two-step) RT-PCR. Components of a replication reaction mixture typically include a nucleic acid template (in this instance the cDNA); a nucleic acid polymerase; and the appropriate nucleotide building blocks needed for nucleic acid synthesis. Nucleic acid replication refers to the polymerization of a nucleic acid whose sequence is determined by, and complementary to, another nucleic acid. DNA replication, as used herein, is synonymous with DNA amplification. Preferably DNA amplification occurs repetitively, thus replicating both strands of the nucleic acid sequence, i.e., DNA complementary to the RNA template, and DNA whose nucleic acid sequence is substantially identical to the RNA template. Repetitive, or cyclic, DNA replication may be advantageously accomplished using a thermostable polymerase in a Polymerase Chain Reaction (PCR).

[0030] PCR is a technique well known in the art. PCR is used to amplify nucleic acids by subjecting a reaction mixture to cycles of: (i) nucleic acid denaturation, (ii) oligonucleotide primer annealizing, and (iii) nucleic acid polymerization. Preferred reaction conditions for amplification comprise thermocycling, i.e., alternating the temperature of the reaction mixture to facilitate each of the steps of the PCR cycle. PCR is typically extended through multiple cycles of denaturation, annealizing and replication, augmented (optionally and preferably) with an initial prolonged denaturation step and a final prolonged extension (polymerization) step. Thermocycling typically occurs within a temperature range of between about 23°C to about 100°C, and preferably between about 37°C to about 95°C. Nucleic acid denaturation typically occurs between about 90°C to about 100°C, preferably about 94°C. Annealizing typically occurs between about 37°C to about 75°C, preferably about 60°C. Polymerization typically occurs between about 55°C to about 80°C, preferably about 72°C. The number of thermocycles varies immensely, depending upon practitioner preference and the quantity of DNA product desired. Preferably, the number of PCR cycles ranges from about 5 to about 99, more preferably greater than about 20 cycles, most preferably about 40 cycles.

Template RNA

[0031] The template RNA may be any ribonucleic acid of interest, known or unknown to the practitioner. Template RNA may be artificially synthesized or isolated from natural sources. Preferably the RNA is mRNA. More preferably the RNA is biologically active or encodes a biologically active polypeptide.

Homopolymeric Nucleic Acids

[0032] The present invention relates to the discovery that homopolymeric nucleic acids serve as an inhibition-relieving agent capable of suppressing, or otherwise reducing, the inhibition of nucleic acid replication by reverse transcriptase, as observed in RT-PCR. A homopolymeric nucleic acid is any oligonucleotide of two or more bases in length composed of a single nucleotide species. Homopolymeric nucleic acids may be ribonucleic acids or deoxyribonucleic acids. Preferred homopolymeric ribonucleic acids include poly(A), poly(C), poly(G) or poly(U). Preferred homopolymeric deoxyribonucleic acids include poly(dA), poly(dC), poly(dG), or poly(dI). The homopolymer used for RT-PCR is an oligonucleotide of two or more bases in length, preferably between 100 to 5000 bases in length, and more preferably between 200 to 2000 bases in length. The amount of nucleic acid homopolymer used for RT-PCR may be greater than 10 ng, preferably between 10 ng to 2000 ng, more preferably between 100 ng to 1000 ng, and most preferably between 200 ng to 600 ng. Preferably, the final concentration of nucleic acid homopolymer in the RT-PCR is between 200 ng/ml to 40000 ng/ml, more preferably between 2000 ng/ml to 20000 ng/ml, and most preferably between 4000 ng/ml to 12000 ng/ml.

Reverse Transcriptases

[0033] Reverse transcriptases useful in the present invention may be any polymerase that exhibits reverse transcriptase activity. Preferred enzymes include those that exhibit reduced RNase H activity. Several reverse transcriptases are known in the art and are commercially available (e.g., from Boehringer Mannheim Corp., Indianapolis, IN; Life Technologies, Inc., Rockville, MD; New England Biolabs, Inc., Beverly, MA; Perkin Elmer Corp., Norwalk, CT; Pharmacia LKB Biotechnology, Inc., Piscataway, NJ; Qiagen, Inc., Valencia, CA; Stratagene, La Jolla, CA). Preferred reverse transcriptases include: Avian Myeloblastosis Virus reverse transcriptase (AMV-RT), Moloney Murine Leukemia Virus reverse transcriptase (MMLV-RT), Human Immunovirus reverse transcriptase (HIV-RT), EIAV-RT, RAV2-RT, *C. hydrogeniformans* DNA Polymerase, *rTth* DNA polymerase, SuperScript I, SuperScript II, and mutants, variants and derivatives thereof. It is to be understood that a variety of reverse transcriptases may be used in the present invention,

including reverse transcriptases not specifically disclosed above, without departing from the scope or preferred embodiments thereof.

DNA Polymerases

- 5 **[0034]** DNA polymerases useful in the present invention may be any polymerase capable of replicating a DNA molecule. Preferred DNA polymerases are thermostable polymerases, which are especially useful in PCR. Thermostable polymerases are isolated from a wide variety of thermophilic bacteria, such as *Thermus aquaticus* (*Taq*), *Thermus brockianus* (*Tbr*), *Thermus flavus* (*Tfl*), *Thermus ruber* (*Tru*), *Thermus thermophilus* (*Tth*), *Thermococcus litoralis* (*Tli*) and other species of the *Thermococcus* genus, *Thermoplasma acidophilum* (*Tac*), *Thermotoga neapolitana* (*Tne*), *Thermotoga maritima* (*Tma*), and other species of the *Thermotoga* genus, *Pyrococcus furiosus* (*Pfu*), *Pyrococcus woesei* (*Pwo*) and other species of the *Pyrococcus* genus, *Bacillus stercophilus* (*Bst*), *Sulfolobus acidocaldarius* (*Sac*) *Sulfolobus solfataricus* (*Sso*), *Pyrodictium occultum* (*Poc*), *Pyrodictium abyssi* (*Pab*), and *Methanobacterium thermoautotrophicum* (*Mth*), and mutants, variants or derivatives thereof.
- 10 **[0035]** Several DNA polymerases are known in the art and are commercially available (e.g., from Boehringer Mannheim Corp., Indianapolis, IN; Life Technologies, Inc., Rockville, MD; New England Biolabs, Inc., Beverly, MA; Perkin Elmer Corp., Norwalk, CT; Pharmacia LKB Biotechnology, Inc., Piscataway, NJ; Qiagen, Inc., Valencia, CA; Stratagene, La Jolla, CA). Preferably the thermostable DNA polymerase is selected from the group of *Taq*, *Tbr*, *Tfl*, *Tru*, *Tth*, *Tli*, *Tac*, *Tne*, *Tma*, *Tih*, *Tfi*, *Pfu*, *Pwo*, *Kod*, *Bst*, *Sac*, *Sso*, *Poc*, *Pab*, *Mth*, *Pho*, ES4, VENT™, DEEPVENT™, and active mutants,
- 15 variants and derivatives thereof. It is to be understood that a variety of DNA polymerases may be used in the present invention, including DNA polymerases not specifically disclosed above, without departing from the scope or preferred embodiments thereof.
- 20

Oligonucleotide Primers

- 25 **[0036]** Oligonucleotide primers useful in the present invention may be any oligonucleotide of two or more nucleotides in length. Preferably, PCR primers are about 15 to about 30 bases in length, and are not palindromic (self-complementary) or complementary to other primers that may be used in the reaction mixture. Primers may be, but are not limited to, random primers, homopolymers, or primers specific to a target RNA template (e.g., a sequence specific primer). Oligonucleotide primers are oligonucleotides used to hybridize to a region of a target nucleic acid to facilitate the polymerization of a complementary nucleic acid. In preferred RT-PCR techniques, primers serve to facilitate reverse transcription of a first nucleic acid molecule complementary to a portion of an RNA template (e.g., a cDNA molecule), and also to facilitate replication of the nucleic acid (e.g., PCR amplification of DNA). Any primer may be synthesized by a practitioner of ordinary skill in the art or may be purchased from any of a number of commercial vendors (e.g., from
- 30 Boehringer Mannheim Corp., Indianapolis, IN; New England Biolabs, Inc., Beverly, MA; Pharmacia LKB Biotechnology, Inc., Piscataway, NJ). It is to be understood that a vast array of primers may be useful in the present invention, including those not specifically disclosed herein, without departing from the scope or preferred embodiments thereof.
- 35

Nucleotide Bases

- 40 **[0037]** Nucleotide bases useful in the present invention may be any nucleotide useful in the polymerization of a nucleic acid. Nucleotides may be naturally occurring, unusual, modified, derivative, or artificial. Nucleotides may be unlabeled, or detectably labeled by methods known in the art (e.g., using radioisotopes, vitamins, fluorescent or chemiluminescent moieties, dioxigenin). Preferably the nucleotides are deoxynucleoside triphosphates, dNTPs (e.g., dATP, dCTP, dGTP, dTTP, dITP, dUTP, α -thio-dNTPs, biotin-dUTP, fluorescein-dUTP, digoxigenin-dUTP, 7-deaza-dGTP). dNTPs are also well known in the art and are commercially available vendors (e.g., from Boehringer Mannheim Corp., Indianapolis, IN; New England Biolabs, Inc., Beverly, MA; Pharmacia LKB Biotechnology, Inc., Piscataway, NJ).
- 45

Buffering Agents and Salt Solutions

- 50 **[0038]** Buffering agents and salts useful in the present invention provide appropriate stable pH and ionic conditions for nucleic acid synthesis, e.g., for reverse transcriptase and DNA polymerase activity. A wide variety of buffers and salt solutions and modified buffers are known in the art that may be useful in the present invention, including agents not specifically disclosed herein. Preferred buffering agents include, but are not limited to, TRIS, TRICINE, BIS-TRICINE, HEPES, MOPS, TES, TAPS, PIPES, CAPS. Preferred salt solutions include, but are not limited to solutions of; potassium acetate, potassium sulfate, ammonium sulfate, ammonium chloride, ammonium acetate, magnesium chloride, magnesium acetate, magnesium sulfate, manganese chloride, manganese acetate, manganese sulfate, sodium chloride, sodium acetate, lithium chloride, and lithium acetate.
- 55

Other Additives Useful in RT-PCR

[0039] Other additives capable of facilitating reverse transcription, replication, and/or a combination of both reactions (e.g., agents for facilitating RT-PCR), other than those disclosed for the first time by this invention, are known in the art. In accordance with the compositions and methods of this invention, one or more of these additives may be incorporated in the present compositions to optimize the generation and replication of nucleic acids from a ribonucleic acid template. Additives may be organic or inorganic compounds. Inhibition-relieving agents useful in the present invention include, but are not limited to, polypeptides such as; human serum albumin, bovine serum albumin (BSA), ovalbumin, albumax, casein, gelatin, collagen, globulin, lysozyme, transferrin, myoglobin, hemoglobin, α -lactalbumin, fumarase, glyceraldehyde-3-phosphate dehydrogenase (GAPDH), amyloglucosidase, carbonic anhydrase, β -lactoglobulin, aprotinin, soybean trypsin inhibitor, trypsinogen, phosphorylase b, myosin, actin, β -galactosidase, catalase, tryptic soy digests, tryptose, lectins, *E. coli* single-stranded binding (SSB) protein, phage T4 gene 32 protein, and the like, or fragments or derivatives thereof. Examples of nonpolypeptide additives include, but are not limited to; tRNA, rRNA, sulfur-containing compounds, acetate-containing compounds, dimethylsulfoxide (DMSO), glycerol, formamide, betain, tetramethylammonium chloride (TMAC), polyethylene glycol (PEG), Tween 20, NP 40, ectoine, and polyols. Preferred additives include DMSO, glycerol, formamide, betain, TMAC, PEG, Tween 20, NP 40, ectoine, polyols, *E. coli* (SSB) protein, Phage T4 gene 32 protein, BSA.

Reagent kits for coupled RT-PCR

[0040] A further aspect of the present invention is a reagent kit for coupled RT-PCR, comprising a reverse transcriptase, a DNA polymerase, and a homopolymeric nucleic acid. In such a reagent kit, said reverse transcriptase, DNA polymerase and homopolymeric nucleic acid may be combined in one container. Alternatively, said reverse transcriptase and DNA polymerase may be combined in one container and said homopolymeric nucleic acid may be in a second container. Alternatively, said reverse transcriptase and said homopolymeric acid may be combined in one container and said DNA polymerase may be in a second container. Alternatively, said DNA polymerase and said homopolymeric nucleic acid may be combined in one container and said reverse transcriptase may be in a second container. Alternatively, all three components may be in separate containers. Such a reagent kit may further comprise a mixture of two or more of dATP, dCTP, dGTP, and dTTP, and/or one or more additives useful in RT-PCR, and/or a buffer. All components mentioned above may be in separate containers, or two or more of them may be combined in one container. Preferably, the two enzymes are combined in one container, either together with the homopolymeric nucleic acid or the homopolymeric nucleic acid being in a separate container, and a buffer solution, a dNTP mix, and one or more useful additive, if present, are provided in the kit in separate containers.

[0041] A further aspect of the present invention is the use of a composition or a reagent kit as above for coupled RT-PCR.

Particular embodiments of the present invention are as follows:

A composition suitable for reverse transcription of RNA, comprising a homopolymeric nucleic acid and at least one member of the group consisting of: a reverse transcriptase, a DNA polymerase, one or more oligonucleotide primers, any one or more nucleotides or derivatives thereof, a buffering agent, a salt, and additives that are useful in RT-PCR.

The composition as above, wherein said homopolymeric nucleic acid is a homopolymeric ribonucleic acid.

The composition as above, wherein said homopolymeric ribonucleic acid is selected from the group consisting of: homopolymeric poly(A), homopolymeric poly(U), homopolymeric poly(C), and homopolymeric poly(G).

The composition as above, wherein said homopolymeric nucleic acid is a homopolymeric deoxyribonucleic acid.

The composition as above, wherein said homopolymeric deoxyribonucleic acid is selected from the group consisting of: homopolymeric poly(dA), homopolymeric poly(dI), homopolymeric poly(dC), and homopolymeric poly(dG).

The composition as above, wherein the composition contains at least 10 ng of homopolymer.

The composition as above, wherein the composition contains between 10 to 2000 ng of homopolymer.

The composition as above, wherein the composition contains between 100 to 1000 ng of homopolymer.

EP 1 050 587 A2

The composition as above, wherein the composition contains between 200 to 600 ng of homopolymer.

The composition as above further comprising a reverse transcriptase.

- 5 The composition as above, wherein said reverse transcriptase is selected from the group consisting of; AMV-RT, M-MLV-RT, HIV-RT, EIAV-RT, RAV2-RT, *C. hydrogenoformans* DNA Polymerase, SuperScript I, SuperScript II, and mutants, variants and derivatives thereof.

The composition as above further comprising a DNA polymerase.

- 10 The composition as above, wherein said DNA polymerase is selected from the group consisting of; *Taq*, *Tbr*, *Tfi*, *Tru*, *Tth*, *Tli*, *Tac*, *Tne*, *Tma*, *Tih*, *Tfi*, *Pfu*, *Pwo*, *Kod*, *Bst*, *Sac*, *Sso*, *Poc*, *Pab*, *Mth*, *Pho*, ES4, VENT™, DEEP-VENT™, and mutants, variants and derivatives thereof.

- 15 The composition as above further comprising one or more oligonucleotide primers.

The composition as above, wherein said primer is selected from the group consisting of; random primers, homopolymers, or primers specific to a target RNA template.

- 20 The composition as above further comprising a buffering agent.

The composition as above, wherein said buffering agent is selected from the group consisting of; TRIS, TRICINE, BIS-TRICINE, HEPES, MOPS, TES, TAPS, PIPES, and CAPS.

- 25 The composition as above further comprising a salt solution.

The composition as above, wherein said salt solution is selected from the group consisting of; potassium chloride, potassium acetate, potassium sulfate, ammonium sulfate, ammonium chloride, ammonium acetate, magnesium chloride, magnesium acetate, magnesium sulfate, manganese chloride, manganese acetate, manganese sulfate, sodium chloride, sodium acetate, lithium chloride, and lithium acetate.

- 30

The composition as above further comprising an additive useful in RT-PCR.

- 35 The composition as above, wherein said additive useful in RT-PCR is selected from the group consisting of; DMSO, glycerol, formamide, betain, tetramethylammonium chloride, PEG, Tween 20, NP 40, ectoine, polyols, *E. coli* SSB protein, Phage T4 gene 32 protein, BSA

The composition as above further comprising one or more nucleotides.

- 40 The composition as above, wherein said nucleotide is a deoxyphosphonucleotide or derivative thereof.

The composition as above, wherein said deoxyphosphonucleotide is selected from the group consisting of; dATP, dCTP, dGTP, dTTP, dITP, dUTP, α -thio-dNTPs, biotin-dUTP, fluorescein-dUTP, and digoxigenin-dUTP.

- 45 A method for generating a nucleic acid from an RNA template comprising the steps:

- a) adding said RNA template to a reaction mixture, said reaction mixture comprising at least one reverse transcriptase and at least one nucleic acid polymerase, or derivatives thereof, and a homopolymeric nucleic acid; and
b) incubating said reaction mixture under conditions sufficient to allow polymerization of a nucleic acid molecule complementary to a portion of said RNA template.

50

A method for generating a nucleic acid from an RNA template comprising the steps:

- a) adding said RNA template to a reaction mixture, said reaction mixture comprising at least one reverse transcriptase, at least one nucleic acid polymerase, one or more ribonucleotide triphosphate or derivatives thereof, and a homopolymeric nucleic acid; and
b) incubating said reaction mixture under conditions sufficient to allow polymerization of a nucleic acid molecule complementary to a portion of said RNA template.

55

[0042] The method as above wherein said nucleic acid is DNA, said polymerase is DNA polymerase, and said ribonucleic acid triphosphate derivative is selected from the group consisting of; dATP, dCTP, dGTP, dTTP, dITP, dUTP, α -thio-dNTPs, biotin-dUTP, fluorescein-dUTP, and digoxigenin-dUTP.

[0043] The method as above, wherein said homopolymeric nucleic acid is a homopolymeric ribonucleic acid.

[0044] The method as above, wherein said homopolymeric ribonucleic acid is selected from the group consisting of; homopolymeric poly(A), homopolymeric poly(U), homopolymeric poly(C), and homopolymeric poly(G).

[0045] The method as above, wherein said homopolymeric nucleic acid is a homopolymeric deoxyribonucleic acid.

[0046] The method as above, wherein said homopolymeric deoxyribonucleic acid is selected from the group consisting of; homopolymeric poly(dA), homopolymeric poly(dI), homopolymeric poly(dC), and homopolymeric poly(dG).

[0047] The method as above, wherein said DNA polymerase is a thermostable polymerase.

[0048] The method as above, wherein said DNA polymerase is selected from the group consisting of; *Taq*, *Tbr*, *Tfi*, *Tru*, *Tth*, *Tli*, *Tac*, *Tne*, *Tma*, *Tih*, *Tfi*, *Pfu*, *Pwo*, *Kod*, *Bst*, *Sac*, *Sso*, *Poc*, *Pab*, *Mth*, *Pho*, ES4, VENT™, DEEPVENT™, and mutants, variants and derivatives thereof.

[0049] The method as above, wherein said reverse transcriptase is selected from the group consisting of; AMV-RT, M-MLV-RT, HIV-RT, EIAV-RT, RAV2-RT, *C. hydrogenoformans* DNA Polymerase, SuperScript I, SuperScript II, and mutants, variants and derivatives thereof.

[0050] The method as above, wherein said reaction mixture further comprises a nucleic acid primer

[0051] The method as above, wherein said primer is complementary to a portion of said RNA template.

[0052] The method as above, wherein said DNA molecule complementary to a portion of said RNA template is amplified by a process of DNA replication.

[0053] The method as above, wherein said process of DNA replication comprises a polymerase chain reaction (PCR).

[0054] The method as above, wherein the combined reactions of DNA generation and amplification occurs in a coupled reaction mixture.

[0055] The method as above, wherein the reaction proceeds between about 23°C to about 100°C.

[0056] The method as above, wherein reverse transcription occurs between about 37°C to about 70°C and said PCR amplification occurs between about 60°C to about 95°C

[0057] The method as above, wherein said PCR amplification comprises at least 25 cycles of denaturation, anneal-ization, and polymerization.

[0058] The method as above, wherein said reverse transcription occurs between about 37°C to about 60°C, dena-turation occurs at about 94°C, anneal-ization occurs at about 60°C, and polymerization occurs at about 72°C.

[0059] It will be readily apparent to one of ordinary skill in the art that other suitable modifications and adaptations of the compositions and methods of the invention described herein are obvious and may be made without departing from the scope of the invention or the disclosed embodiments thereof. Having now described the present invention in detail, the same will be more clearly understood by reference to the following examples, which are included herewith for purposes of illustration only and are not intended to be limiting of the invention.

EXAMPLE 1: Inhibition of Nucleic Acid Amplification by Reverse Transcriptase

[0060] To study the inhibiting effect of reverse transcriptases on nucleic acid replication, 10^3 and 10^2 copies of a plasmid containing a 500 bp fragment of the RanBP2 gene were used as template for PCR amplification. The reactions were performed in a final volume of 50 μ l. The reaction mixture comprised 50 mM KCl, 10 mM Tris-HCl, pH 8.3, 0.001% gelatine, 1.5 mM $MgCl_2$, 2 mM DTT, 200 μ M of each dNTP, 0.16 μ M sense and antisense primers (SEQ ID NOs: 1 and 2, respectively), 2.5 Units HotStarTaq™ DNA polymerase (QIAGEN), 10 Units of recombinant RNasin (Promega™) and increasing amounts Omniscript RT (QIAGEN) from 0 to 5 units. To simulate a coupled RT-PCR, the PCR was preceded by an incubation step for 30 minutes at 37°C ("pseudo RT"). The PCR program consisted of an initial RT denatura-tion/HotStarTaq activation step (15 min, 95°C), 45 cycles (30 sec, 94°C; 1 min, 60°C, 1 min, 72°C) and a final extension step (10 min, 72°C). PCR was performed in a Biometra Uno II thermocycler.

[0061] The samples (in triplicate) were analyzed on a 1% agarose gel, and product yield was densitometrically determined. Product yield in the samples without reverse transcriptase was set to 100%.

Table 1

Inhibitory effect of Omniscript RT on PCR							
		Omniscript RT (Units)					
		0	0.2	0.5	1.0	2.0	5.0
Template (copy number)	10 ³	100%	84%	47%	25%	14%	2%
	10 ²	100%	90%	32%	20%	8%	1%

[0062] The results clearly demonstrate that the presence of even little amounts of reverse transcriptase dramatically reduced DNA replication efficiency when only a limiting amount of template DNA was present. Comparable experiments using genomic DNA as template confirmed the inhibiting effect of Omniscript RT on PCR, especially when minute amounts of template DNA were used.

EXAMPLE 2: Reduction of Reverse Transcriptase inhibition of Nucleic Acid Amplification by Homopolymeric poly(A)

[0063] To demonstrate the discovery that the inhibitory effect of reverse transcriptase on DNA replication is relieved by the use of homopolymeric poly(A), increasing amounts of poly(A) (commercially obtained from Pharmacia) were added to samples containing various amounts of total RNA from HeLa cells. A 626 bp fragment of the β -actin mRNA from HeLa total RNA was selected as the RNA template for RT-PCR.

[0064] Reactions were conducted in a final volume of 50 μ l. Two reaction mixtures comprised 50mM KCl, 10mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM MgCl₂, 2 mM DTT, 400 μ M each dNTP, 0.4 μ M sense and antisense primers (SEQ ID NOs: 3 and 4, respectively), 5 Units HotStar Taq DNA polymerase (QIAGEN), 0.5 Units Omniscript RT and 1 and 10 ng of HeLa total RNA, respectively. Reactions were performed in the presence of 0 ng, 100 ng, 200 ng and 400 ng of poly(A) homopolymer. Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation/HotStarTaq activation step (15 min, 95°C), 40 cycles (1 min, 94°C; 1 min, 60°C, 1 min, 72°C) and a final extension step (10 min 72°C). RT-PCR was performed in an Eppendorf gradient cyclor.

[0065] The samples (in duplicate) were analyzed on a 1% agarose gel. A clear improvement in efficiency and specificity upon addition of poly(A) was observed for both the 1 and 10 ng of total HeLa RNA samples. The yield of RT-PCR significantly increased with increasing amounts of poly(A). Figure 1 shows the effect of poly(A) when using 1 ng of HeLa total RNA as template. The same effect was observed when 10 ng of total HeLa RNA were used as template (data not shown).

The results demonstrate that the addition of poly(A) homopolymer reduces the inhibiting effect of Omniscript RT (QIAGEN) on HotStar Taq DNA Polymerase (QIAGEN), thereby contributing to higher sensitivity and specificity in the one-step RT-PCR.

EXAMPLE 3: Reduction of Reverse Transcriptase inhibition of Nucleic Acid Amplification by Other Homopolymeric Ribonucleic Acids

[0066] The experiments described in Example 2 were replicated, using poly(U) and poly(C) in place of poly(A), to demonstrate that the reduction of the inhibitory effect of reverse transcriptase on DNA replication, due to the presence of a homopolymeric ribonucleic acid, is a characteristic common to homopolymeric ribonucleic acids generally.

[0067] The same 626 bp fragment of β -actin mRNA was selected as the template RNA. Reactions were conducted in a final volume of 50 μ l. The reaction mixture comprised 50 mM KCl, 10 mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM MgCl₂, 2 mM DTT, 400 μ M each dNTP, 0.4 μ M sense and antisense primers (SEQ ID NOs: 3 and 4, respectively), 5 Units HotStar Taq DNA polymerase (QIAGEN), 0.5 Units Omniscript RT (QIAGEN) and 1 ng of HeLa total RNA. Reactions were performed in the presence of 0 ng, 100 ng, 200 ng and 400 ng of poly(U) (Pharmacia) and poly(C) (Pharmacia) homopolymer, respectively. Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation/HotStarTaq activation step (15 min, 95°C), 40 cycles (1 min, 94°C; 1 min, 60°C, 1 min, 72°C) and a final extension step (10 min, 72°C). RT-PCR was performed in an Eppendorf gradient cyclor.

[0068] The samples (in triplicate) were analyzed on a 1% agarose gel. Figures 2a and 2b show the effect of poly(U) and poly(C), respectively, when using 1 ng of HeLa total RNA as template and 0 ng, 100 ng and 400 ng of the respective homopolymeric ribonucleic acid. The results of this series of experiments confirm that the addition of homopolymeric ribonucleic acids suppress or reduce the inhibitory effect of reverse transcriptase on DNA replication.

EXAMPLE 4: Reduction of Reverse Transcriptase inhibition of Nucleic Acid Amplification by Homopolymeric Deoxyribonucleic Acids

[0069] The experiments, as described in Examples 2 and 3, were repeated using a homopolymeric deoxyribonucleic acid (poly(dA)), to demonstrate that the reduction of the inhibitory effect of reverse transcriptase on DNA replication, is a feature common to homopolymeric nucleic acids generally.

[0070] The same 626 bp fragment of β -actin mRNA was selected as the template RNA. Reactions were conducted in a final volume of 50 μ l. The reaction mixture comprised 50 mM KCl, 10mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM $MgCl_2$, 2 mM DTT, 400 μ M each dNTP, 0.4 μ M sense and antisense primers (SEQ ID NOs: 3 and 4 respectively), 5 Units HotStar Taq DNA polymerase (QIAGEN), 0.5 Units Omniscript RT (QIAGEN) and 1 ng of HeLa total RNA. Reactions were performed in the presence of 0 ng, 100 ng, 200 ng and 400 ng of poly(dA) homopolymer (Pharmacia). Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation/HotStar Taq activation step (15 min, 95°C), 40 cycles (1 min, 94°C; 1 min, 60°C, 1 min, 72°C) and a final extension step (10 min, 72°C). RT-PCR was performed in an Eppendorf gradient cyclor.

[0071] The samples (in triplicate) were analyzed on a 1% agarose gel. Figure 3 illustrates the addition of poly(dA) substantially improves RT-PCR. These results confirm that the inhibitory effect of reverse transcriptase on DNA replication is reduced by addition of homopolymeric nucleic acids, generally (deoxyribonucleic acids as well as ribonucleic acids).

EXAMPLE 5: Reduction of Reverse Transcriptase inhibition of Nucleic Acid Amplification by Homopolymeric Nucleic Acids Under Varying Reverse Transcriptase Conditions

[0072] To demonstrate that the effectiveness of the present invention to reduce the inhibitory effect of reverse transcriptase on DNA replication is independent of the reverse transcriptase used, experiments of the earlier examples were repeated using a variety of enzymes that exhibit reverse transcriptase activity.

[0073] The same 626 bp fragment of β -actin mRNA was selected as the template RNA. Reactions were conducted in a final volume of 50 μ l. The reaction mixture comprised 50 mM KCl, 10 mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM $MgCl_2$, 2 mM DTT, 400 μ M of each dNTP, 0.4 μ M sense and antisense primers (SEQ ID NOs: 3 and 4 respectively), 5 Units HotStar Taq DNA polymerase (QIAGEN), and 1 and 10 ng of HeLa total RNA, respectively. Three reverse transcriptases were tested. The reaction mixture contained either; 10 Units AMV-RT (Roche Molecular Biochemicals), 50 Units M-MLV-RT (Promega, recombinant), or 50 Units M-MLV, RNaseH⁻-RT (Superscript II RT, LTI, recombinant). Reactions were performed in the presence of 0 ng, 100 ng, 200 ng and 400 ng of poly(A) homopolymer. Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation/HotStarTaq activation step (15 min, 95°C), 40 cycles (1 min, 94°C; 1 min, 60°C, 1 min, 72°C) and a final extension step (10 min, 72°C). RT-PCR was performed in a Biometra Uno II thermal cyclor.

[0074] A clear increase in sensitivity (i.e., reduced inhibition of PCR in the one-step RT-PCR) was observed for all three reverse transcriptases tested. Figure 4 summarizes the results (in triplicate) with β -actin as the template in the absence (-) or presence of 400 ng poly(A) homopolymer (+). For AMV-RT and MMLV-RT, results with 1 ng of total RNA as template are shown. For M-MLV, RNaseH⁻ RT, the results obtained with 10 ng template RNA are shown.

[0075] These results confirm the effectiveness of the present invention to reduce the inhibitory effect of reverse transcriptase on DNA replication, independent of the reverse transcriptase used in the reaction mixture. The present invention can be used as a general optimization component of RT-PCR, irrespective of the reverse transcriptase used.

EXAMPLE 6: Reduction of Reverse Transcriptase Inhibition of Nucleic Acid Amplification by Homopolymeric Nucleic Acids Under Varying DNA Polymerase, Buffer and Salt Conditions

[0076] To further demonstrate the general utility of the present invention to reduce the inhibitory effect of reverse transcriptase on DNA replication in a variety of reaction mixtures, trials were conducted employing different DNA polymerases, buffering agents and salt solutions.

[0077] Two DNA polymerases were tested; unmodified, recombinant Taq DNA-Polymerase (AmpliTaq, Perkin Elmer), and recombinant Pfu DNA polymerase (Stratagene). One trial series employed 5 Units AmpliTaq DNA Polymerase in a reaction mixture comprising 50 mM KCl, 10 mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM $MgCl_2$, 400 μ M of each dNTP, and 1 ng of HeLa total RNA. The second trial series employed 5 units Pfu DNA polymerase in a reaction mixture comprising 20 mM Tris-HCl, pH 8.8, 10 mM KCl, 10 mM $(NH_4)_2SO_4$, 2 mM $MgSO_4$, 0.1% Triton X-100, 0.1 mg/ml BSA, 200 μ M of each dNTP, and 1 ng or 100 ng of HeLa total RNA. All reaction mixtures contained 2 mM DTT, 0.4 μ M sense and antisense primers (SEQ ID NOs: 3 and 4, respectively), and 0.5 units Omniscript RT. Reactions were performed in a final volume of 50 μ l.

[0078] Reactions were performed in the presence of 0 ng, 100 ng, 200 ng and 400 ng of poly(A) homopolymer.

Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation step (5 min, 94°C), 40 cycles (1 min, 94°C; 1 min, 60°C, 1 min, 72°C) and a final extension step (10 min, 72°C). RT-PCR was performed in a Biometra Uno II thermal cycler.

[0079] Results confirm improved nucleic acid amplification occurred in the presence of the homopolymeric nucleic acid, regardless of DNA polymerase, buffering agent or salt solution used in the final reaction mixture. Figures 5a and 5b depict the results using 1 ng of total HeLa RNA in the absence (-), or presence (+) of 400 ng poly(A). The inhibition-relieving effect of homopolymeric nucleic acids in coupled RT-PCR is not limited to a particular thermostable DNA polymerase, a particular buffering agent or salt solution. It has, therefore, been demonstrated that the present invention is useful under varying ionic conditions, different pH values, and can be applied to various DNA polymerases.

EXAMPLE 7: Reduction of Reverse Transcriptase inhibition of Nucleic Acid Amplification by Homopolymeric Nucleic Acids On Different Target Ribonucleic Acids

[0080] To demonstrate that the effectiveness of the present invention to reduce the inhibitory effect of reverse transcriptase on DNA replication is independent of the target RNA sought and the oligonucleotide primers used, trials were conducted seeking to reverse transcribe and replicate a different gene fragment. For these trials a 690 bp fragment of the α -catenin mRNA from total HeLa RNA was selected as the RNA template for RT-PCR.

[0081] Reactions were conducted in a final volume of 50 μ l. The reaction mixture comprised 50mM KCl, 10 mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM MgCl₂, 2 mM DTT, 400 μ M of each dNTP, 0.4 μ M sense and antisense primers (SEQ ID NOs: 5 and 6, respectively), 5 Units HotStar Taq DNA polymerase (QIAGEN), 0.5 Units Omniscript RT and 5 ng of total HeLa RNA. Reactions were performed in the presence of either 0 ng, or 500 ng of poly(A) homopolymer. Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation/HotStar Taq activation step (15 min, 95°C), 40 cycles (1 min, 94°C; 1 min, 50°C, 1 min, 72°C) and a final extension step (10 min, 72°C). RT-PCR was performed in an MJR PTC 200 thermal cycler.

[0082] The samples (in triplicate) were analyzed on a 1% agarose gel. A clear improvement in efficiency and specificity upon addition of poly(A) was observed for reverse transcription and amplification of the α -catenin template from 5 ng of HeLa total RNA (see Fig 6). The results confirm the inhibition-relieving effect of homopolymeric nucleic acids in coupled RT-PCR is specific to a particular template RNA or use of a particular oligonucleotide primer.

[0083] The results demonstrate that the addition of poly(A) homopolymer reduces the inhibiting effect of Omniscript RT (QIAGEN) on HotStarTaq DNA Polymerase (QIAGEN), thereby contributing to higher sensitivity and specificity in the one-step RT-PCR.

EXAMPLE 8: Reduction of Reverse Transcriptase inhibition of Nucleic Acid Amplification by Homopolymeric Nucleic Acids Compared to Other Nucleic Acids

[0084] The previous examples have shown, conclusively, the positive effect of homopolymeric nucleic acids in reducing the inhibitory effect of reverse transcriptase on DNA replication independent of the reverse transcriptase, DNA polymerase, oligonucleotide primer, RNA template, buffering agent, or salt solution used in the reaction mixture.

[0085] To demonstrate the unique benefits of the present invention to reduce the inhibitory effect of reverse transcriptase on DNA replication, simultaneous trials were conducted comparing the effect on homopolymeric nucleic acids to other nucleic acids, specifically non-homologous RNA molecules.

[0086] The 690 bp fragment of the α -catenin mRNA was selected as the RNA template. Reactions were conducted in a final volume of 50 μ l. The reaction mixture comprised 50 mM KCl, 10 mM Tris-HCl, pH 8.3, 0.001% gelatine, 2.5 mM MgCl₂, 2 mM DTT, 400 μ M of each dNTP, 0.4 μ M sense and antisense primers (SEQ ID NOs: 5 and 6, respectively), 5 Units HotStar Taq DNA polymerase (QIAGEN), 0.5 Units Omniscript RT (QIAGEN), and 5 ng of HeLa total RNA.

[0087] The trial series consisted of reaction mixtures containing either no additional nucleic acid, 500 ng of poly(A) (Pharmacia), 500 ng tRNA (Sigma), or 500 ng *E. coli* 16S, 23 S rRNA (Roche Molecular Biochemicals). Reverse transcription was done at 42°C for 30 min. The PCR program consisted of an initial RT denaturation/HotStar Taq activation step (15 min, 95°C), 40 cycles (1 min, 94°C; 1 min, 50°C, 1 min, 72°C) and a final extension step (10 min, 72°C). RT-PCR was performed in an MJR PTC 200 thermal cycler.

[0088] Figure 6 shows the comparative results of RT-PCR in the absence (-) or presence (+) of additional nucleic acids. When, except for template RNA, no further RNA was present in the reaction mixture, hardly any DNA product could be detected. Similar negative results were obtained when 500 ng tRNA was added to the reaction mixture, thereby corroborating the observations of Chandler *et al.* (1998) who could not confirm the positive effect of tRNA in one-step RT-PCR postulated by Sellner *et al.* (1992). The addition of 500 ng rRNA to the reaction led to the amplification of the desired nucleic acid product. A second unspecific band in the size range of 400 bp also appeared, however. This result clearly illustrates the clear disadvantages of using rRNA as inhibition-relieving agent based on theoretical

considerations that such molecules serve as potential targets for mispriming. Only the utilization of the homopolymeric nucleic acid resulted in a clear improvement of RT-PCR, confirming the superiority of the compositions and methods of the present invention to reduce the inhibitory effect of reverse transcriptase on DNA replication.

[0089] The present invention incorporates by reference in their entirety techniques well known in the field of molecular biology. These techniques include, but are not limited to, techniques described in the following publications:

- Old, R.W. & S.B. Primrose, Principles of Gene Manipulation: An Introduction To Genetic Engineering (3d Ed. 1985) Blackwell Scientific Publications, Boston. Studies in Microbiology; V.2:409 pp. (ISBN 0-632-01318-4).
 Sambrook, J. et al. eds., Molecular Cloning: A Laboratory Manual (2d Ed. 1989) Cold Spring Harbor Laboratory Press, NY. Vols. 1-3. (ISBN 0-87969-309-6).
 Winnacker, E.L. From Genes To Clones: Introduction To Gene Technology (1987) VCH Publishers, NY (translated by Horst Ibelgaufs). 634 pp. (ISBN 0-89573-614-4).

References

[0090]

- Chadwick et al. Comparison of three RNA amplification methods as sources of DNA for sequencing. BioTechniques 25(5): 818-822 (Nov. 1998).
 Chandler et al. Reverse transcriptase (RT) inhibition of PCR at low concentrations of template and its implications for quantitative RT-PCR. Appl. and Environm. Microbiol. 64(2):669-677 (Feb. 1998).
 Compton, J. Nucleic acid sequence-based amplification. Nature 350:91-92 (1991).
 Cusi et al. Comparison of M-MLV reverse transcriptase and *Tth* polymerase activity in RT-PCR of samples with low virus burden. BioTechniques 17:1034-1036 (1994).
 Guatelli et al. Isothermal, in vitro amplification of nucleic acids by a multienzyme reaction modeled after retroviral replication. PNAS USA 87:1874-1878 (1990).
 Juhasz et al. Sensitivity of tyrosinase mRNA detection by RT-PCR: *rth* DNA polymerase vs. MMLV-RT and Ampli-
 itaq® polymerase. BioTechniques 20:592-600 (1996).
 Mallet et al. Continuous RT-PCR using AMV-RT and *Taq* DNA polymerase: characterization and comparison to
 uncoupled procedures. BioTechniques 18:678-687 (1995).
 Myers and Gelfand. Reverse transcription and DNA amplification by a *Thermus thermophilus* DNA polymerase.
Biochemistry 30:7661-7666 (1991).
 Saiki et al. Enzymatic amplification of beta-globin genomic sequences and restriction site analysis for diagnosis of
 sickle cell anaemia. Science 230:1350-1354 (1985).
 Sellner and Turbett. Comparison of three RT-PCR methods. BioTechniques 25(2):230-234 (Aug. 1998).
 Sellner et al. Reverse transcriptase inhibits *Taq* polymerase activity. Nucl. Acids Res. 20(7):1487-1490 (Apr. 11,
 1992).
 Sellner, et al. A one-tube, one manipulation RT-PCR reaction for detection of Ross River virus. J. Virol. Methods
 40:255-264 (1992).
 Sooknanan et al. Fidelity of nucleic acid amplification with avian myeloblastosis virus reverse transcriptase and T7
 RNA polymerase. BioTechniques 17:1077-1085 (1994).
 Xu, et al. The application of polymerase chain reaction to the detection of rotaviruses in faeces. J. Virol. Methods
 27:29-38 (1990).
 Young et al. Detection of hepatitis C virus RNA by a combined reverse transcription-polymerase chain reaction
 assay. J. Clin. Microbiol. 31:882-886 (1993).

[0091] Each of the publications mentioned herein is incorporated by reference.

Sequences referred to in the examples

[0092]

SEQ ID NO.: 1: ATGTTAGTGA AGAAGAGGAG GATG
 SEQ ID NO.: 2: CTTGGCTTGT AGAATCTGTA TCAA
 SEQ ID NO.: 3: CCTTGCCTTT GCCGATCC
 SEQ ID NO.: 4: GGATCTTCAT GAGGTAGTCA GTC
 SEQ ID NO.: 5: AGCTGAAAGT TGTGGAAGAT G
 SEQ ID NO.: 6: GGAGCTGTCT ACGCAAGTC

Sequence Listing

5 SEQUENCE LISTING

<110> Missel, Andreas
 Loeffert, Dirk
 Kang, Jie

10 <120> The Generation and Amplification of Nucleic Acids
 from Ribonucleic Acids

<130> QGN-007.0 US

<140>

<141>

<150>

<151>

15 <160> 6

<170> Microsoft Word 97

<210> 1

20 <211> 24

<212> DNA

<213> Artificial sequence

<220>

<221> Primer_bind

<222>

25 <223> 5' primer for RANBP2

<400> 1

atgttagtga agaagaggag gatg
 24

30 <210> 2

<211> 24

<212> DNA

<213> Artificial sequence

<220>

35 <221> Primer_bind

<222>

<223> 3' primer for RANBP2

<400> 2

40 cttggcttgt agaatctgta tcaa
 24

<210> 3

45 <211> 18

<212> DNA

<213> Artificial sequence

<220>

<221> Primer_bind

<222>

50 <223> 5'primer for beta-actin

<400> 3

55

Sequence Listing

ccttgccttt gccgatcc
18

5

<210> 4
<211> 23
<212> DNA
10 <213> Artificial sequence
<220>
<221> Primer_bind
<222>
15 <223> 3'primer for beta-actin
<400> 4

ggatcttcat gaggtagtca gtc
23

20

<210> 5
<211> 21
<212> DNA
25 <213> Artificial sequence
<220>
<221> Primer_bind
<222>
30 <223> 5'primer for alpha-catenin
<400> 5

agctgaaagt tgtggaagat g
21

35

<210> 6
<211> 19
<212> DNA
40 <213> Artificial sequence
<220>
<221> Primer_bind
<222>
45 <223> 3'primer for alpha-catenin
<400> 6

ggagctgtct acgcaagtc
19

50

Claims

- 55 1. A composition suitable for reverse transcription of RNA, comprising a homopolymeric nucleic acid and at least one member of the group consisting of: a reverse transcriptase, a DNA polymerase, one or more oligonucleotide primers, any one or more nucleotides or derivatives thereof, a buffering agent, a salt, and additives that are useful in RT-PCR.

EP 1 050 587 A2

2. The composition of claim 1, wherein said homopolymeric nucleic acid is a homopolymeric ribonucleic acid, or a homopolymeric deoxyribonucleic acid, or a mixture of both.
3. The composition of claim 2, wherein said homopolymeric ribonucleic, deoxyribonucleic acid, or mixture thereof is homopolymeric poly(A), homopolymeric poly(U), homopolymeric poly(C), homopolymeric poly(G), homopolymeric poly(dA), homopolymeric poly(dI), homopolymeric poly(dC), or homopolymeric poly(dG), or a mixture of two or more of these.
4. The composition of any one of claims 1 to 3, wherein the composition contains at least 10 ng or at least 200 ng/ml of homopolymer.
5. The composition of claim 4, wherein the composition contains between 10 to 2000 ng or between 200 and 40000 ng/ml of homopolymer.
6. The composition of any one of claims 1 to 5 comprising a reverse transcriptase.
7. The composition of claim 6, wherein said reverse transcriptase is AMV-RT, M-MLV-RT, HIV-RT, EIAV-RT, RAV2-RT, *C. hydrogeniformans* DNA Polymerase, SuperScript I, SuperScript II, and/or mutants, variants and derivatives thereof.
8. The composition of any one of claims 1 to 7 comprising a DNA polymerase.
9. The composition of claim 8, wherein said DNA polymerase is *Taq*, *Tbr*, *Tfi*, *Tru*, *Tth*, *Tli*, *Tac*, *Tne*, *Tma*, *Tih*, *Tfi*, *Pfu*, *Pwo*, *Kod*, *Bst*, *Sac*, *Sso*, *Poc*, *Pab*, *Mth*, *Pho*, ES4, VENT™, DEEPVENT™ DNA polymerase, and/or mutants, variants and derivatives thereof.
10. The composition of any one of claims 1 to 9 comprising a reverse transcriptase and a DNA polymerase.
11. The composition of any one of claims 1 to 10 further comprising one or more oligonucleotide primers, said primers preferably being random primers, homopolymers, and/or primers specific to a target RNA template.
12. The composition of any one of claims 1 to 11 further comprising a buffering agent, said buffering agent preferably being TRIS, TRICINE, BIS-TRICINE, HEPES, MOPS, TES, TAPS, PIPES, and/or CAPS..
13. The composition of any one of claims 1 to 12 further comprising a salt solution, said salt solution preferably being potassium chloride, potassium acetate, potassium sulfate, ammonium sulfate, ammonium chloride, ammonium acetate, magnesium chloride, magnesium acetate, magnesium sulfate, manganese chloride, manganese acetate, manganese sulfate, sodium chloride, sodium acetate, lithium chloride, and/or lithium acetate.
14. The composition of any one of claims 1 to 13 further comprising an additive useful in RT-PCR, said additive useful in RT-PCR preferably being DMSO, glycerol, formamide, betain, tetramethylammonium chloride, PEG, Tween 20, NP 40, ectoine, polyols, *E. coli* SSB protein, Phage T4 gene 32 protein, and/or BSA.
15. The composition of any one of claims 1 to 14 further comprising one or more nucleotides, said nucleotides preferably being dATP, dCTP, dGTP, dTTP, dITP, dUTP, α -thio-dNTPs, biotin-dUTP, fluorescein-dUTP, and/or digoxigenin-dUTP.
16. A method for generating a nucleic acid from an RNA template comprising the steps:
 - a) adding said RNA template to a reaction mixture, said reaction mixture comprising at least one reverse transcriptase and at least one nucleic acid polymerase, or derivatives thereof, and a homopolymeric nucleic acid; and
 - b) incubating said reaction mixture under conditions sufficient to allow polymerization of a nucleic acid molecule complementary to a portion of said RNA template.
17. A method for generating a nucleic acid from an RNA template comprising the steps:
 - a) adding said RNA template to a reaction mixture, said reaction mixture comprising at least one reverse tran-

scriptase, at least one nucleic acid polymerase, one or more ribonucleotide triphosphate or derivatives thereof, and a homopolymeric nucleic acid; and

b) incubating said reaction mixture under conditions sufficient to allow polymerization of a nucleic acid molecule complementary to a portion of said RNA template.

- 5 18. The method of claim 17 wherein said nucleic acid is DNA, said polymerase is DNA polymerase, and said ribonucleic acid triphosphate derivative is dATP, dCTP, dGTP, dTTP, dITP, dUTP, α -thio-dNTPs, biotin-dUTP, fluorescein-dUTP, and/or digoxigenin-dUTP.
- 10 19. The method of any one of claims 16 to 18 wherein said homopolymeric nucleic acid is a homopolymeric ribonucleic acid, or a homopolymeric deoxyribonucleic acid, or a mixture of both.
20. The method of claim 19, wherein said homopolymeric nucleic acid is homopolymeric poly(A), homopolymeric poly(U), homopolymeric poly(C), homopolymeric poly(G), homopolymeric poly(dA), homopolymeric poly(dI), homopolymeric poly(dC), and/or homopolymeric poly(dG).
- 15 21. The method of any one of claims 16 to 20, wherein said DNA polymerase is a thermostable polymerase, said thermostable polymerase preferably being *Taq*, *Tbr*, *Tfi*, *Tru*, *Tth*, *Tli*, *Tac*, *Tne*, *Tma*, *Tih*, *Tfi*, *Pfu*, *Pwo*, *Kod*, *Bst*, *Sac*, *Sso*, *Poc*, *Pab*, *Mth*, *Pho*, ES4, VENT™, DEEPVENT™ DNA polymerase, and/or mutants, variants and derivatives thereof.
- 20 22. The method of any one of claims 16 to 21, wherein said reverse transcriptase is AMV-RT, M-MLV-RT, HIV-RT, EIAV-RT, RAV2-RT, *C. hydrogenoformans* DNA Polymerase, SuperScript I, SuperScript II, and/or mutants, variants and derivatives thereof.
- 25 23. The method of any one of claims 16 to 22, wherein said reaction mixture further comprises a nucleic acid primer, said primer preferably being complementary to a portion of said RNA template.
24. The method of claim 23, wherein said DNA molecule complementary to a portion of said RNA template is amplified by a process of DNA replication, said process of DNA replication preferably comprising a polymerase chain reaction (PCR).
- 30 25. The method of claim 24, wherein the combined reactions of DNA generation and amplification occurs in a coupled reaction mixture.
- 35 26. The method of claim 25, wherein the reaction proceeds between about 23°C to about 100°C.
27. The method of claim 26, wherein reverse transcription occurs between about 37°C to about 70°C and said PCR amplification occurs between about 60°C to about 95°C.
- 40 28. The method of any one of claims 24 to 27, wherein said PCR amplification comprises at least 25 cycles of denaturation, annealization, and polymerization.
29. The method of any one of claim 28, wherein said reverse transcription occurs between about 37°C to about 60°C, denaturation occurs at about 94°C, annealization occurs at about 60°C, and polymerization occurs at about 72°C.
- 45 30. Reagent kit for coupled RT-PCR, comprising at least a reverse transcriptase, a DNA polymerase, and a homopolymeric nucleic acid.
- 50 31. Reagent kit of claim 30, wherein said reverse transcriptase, DNA polymerase and homopolymeric nucleic acid are combined in one container, or wherein said reverse transcriptase and DNA polymerase are combined in one container and said homopolymeric nucleic acid is in a second container.
- 55 32. Reagent kit of claim 30 or 31, further comprising a mixture of two or more of dATP, dCTP, dGTP, and dTTP, and/or one or more additives useful in RT-PCR, and/or a buffering agent.
33. Use of a composition according to any one of claims 1 to 15 or a reagent kit according to any one of claims 30 to 32 for coupled RT-PCR.

poly(A) homopolymer

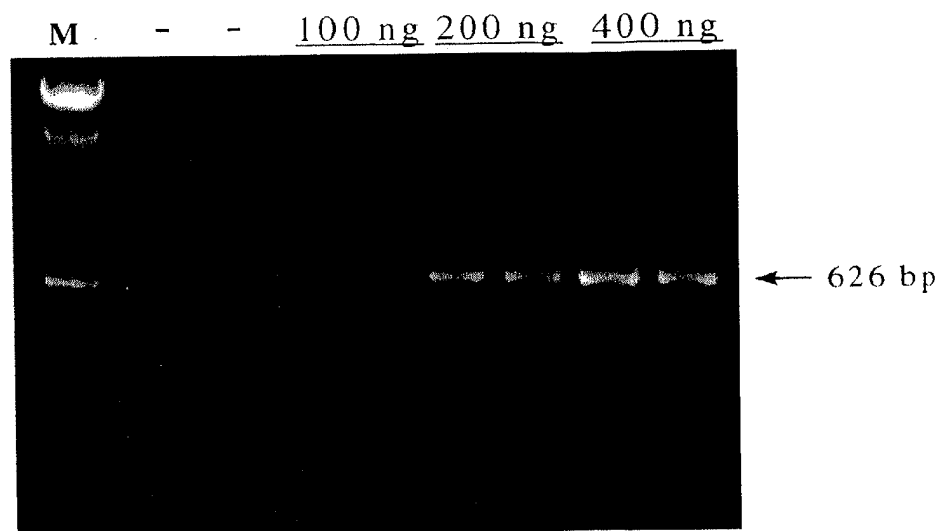


Figure 1: Relief of the inhibiting effect of reverse transcriptase on PCR by homopolymeric poly(A).

FIGURE 1

poly(U) homopolymer

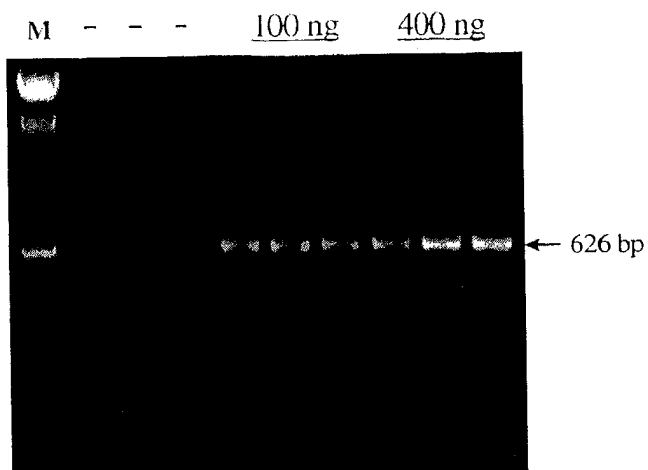


Figure 2A: Relief of the inhibiting effect of reverse transcriptase on PCR by homopolymeric poly(U).

poly(C) homopolymer

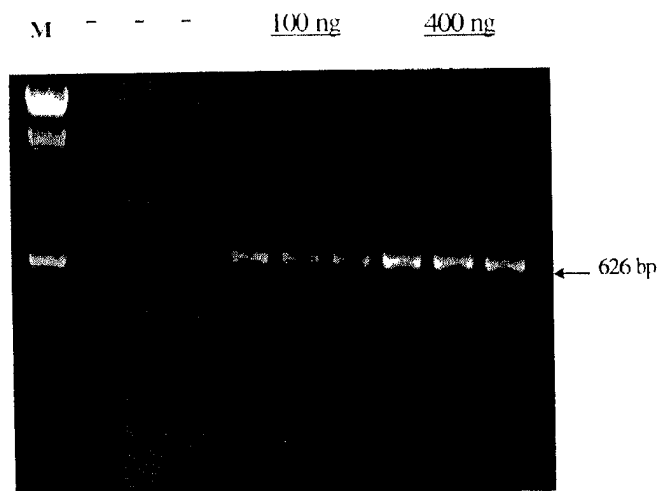


Figure 2B: Relief of the inhibiting effect of reverse transcriptase on PCR by homopolymeric poly(C).

FIGURE 2

EP 1 050 587 A2

poly(dA) homopolymer

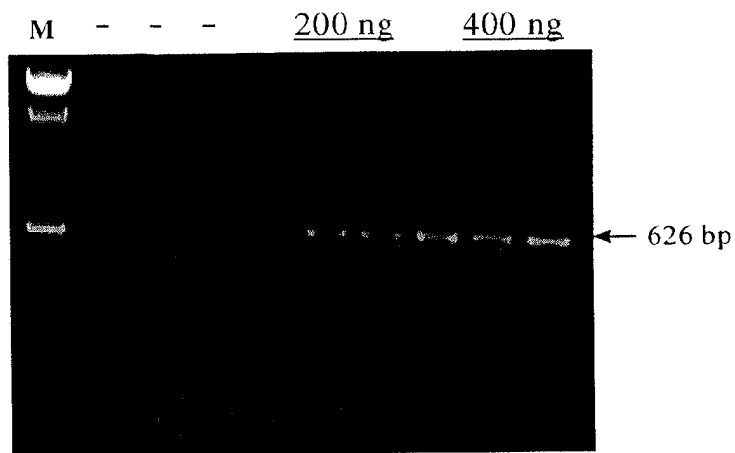


Figure 3: Relief of the inhibiting effect of reverse transcriptase on PCR by homopolymeric poly(dA).

FIGURE 3

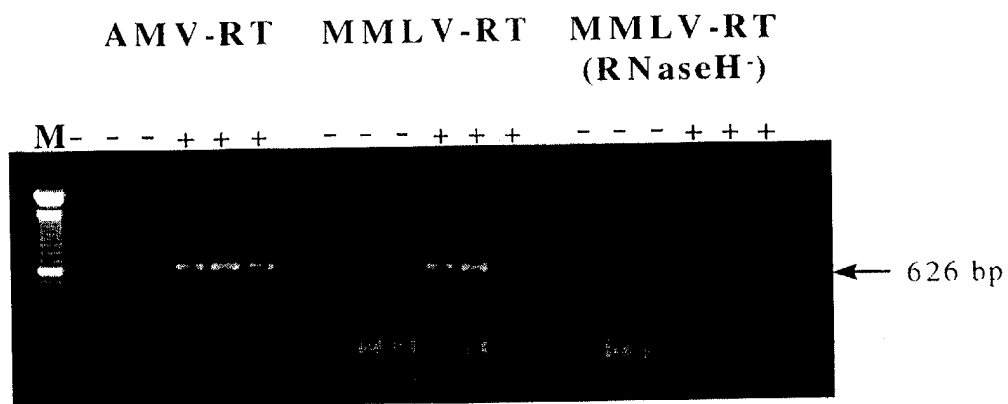


Figure 4: Relief of the inhibiting effect on PCR by poly(A) homopolymer. Effect on different reverse transcriptases.

FIGURE 4

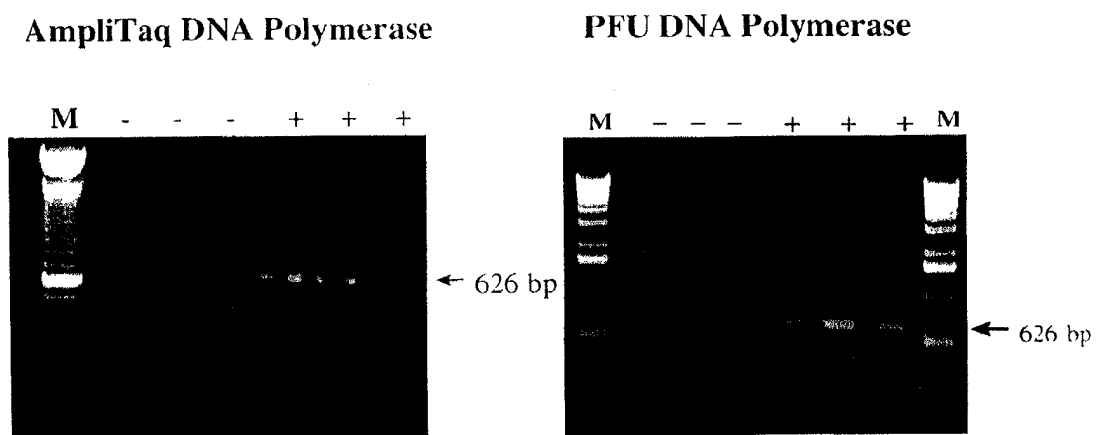


Figure 5: Relief of the inhibiting effect on PCR by poly(A) homopolymer. Effect on different thermostable DNA polymerases and under different buffer conditions.

5a

5b

FIGURE 5

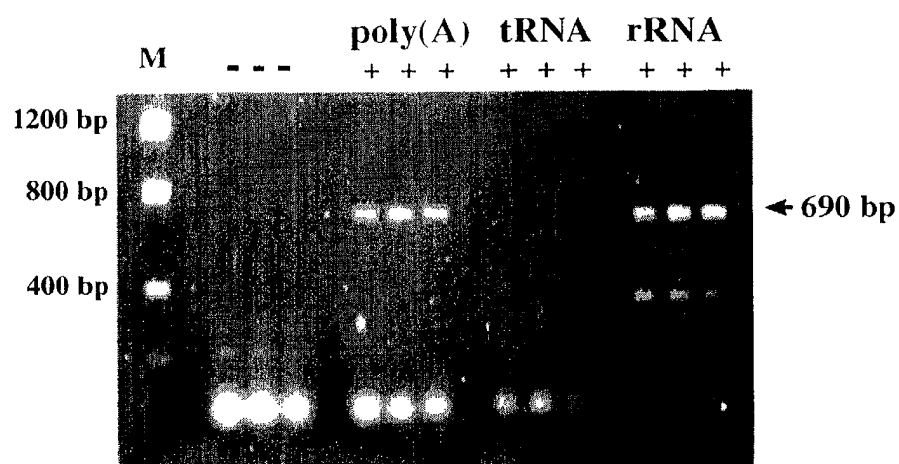


Figure 6: Comparison of different nucleic acids on RT-PCR inhibition relief

FIGURE 6